




N-95 Particulate Respirators

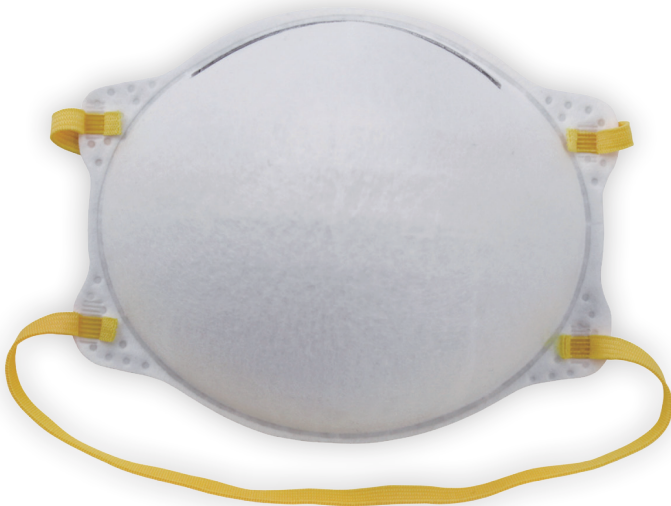

Item No. PRN95240*

Materials

	PRN95240*
Head Straps	•
Nose Foam	•
Nose Clip	•
Filter Media	•

*This model is also cleared by the Food and Drug Administration for use as a Surgical N95 Respirator, FDA 510(k): K020474.

- Materials are latex-free and do not contain any nature rubber latex.

Features and Benefits

Our particulate respirators are designed to protect workers from exposing to solid substances free of oil. The cone style of N95 mask is the most widely used mask in all work environments. These models have been designed to easily meet NIOSH's stringent approval requirements for breathing resistance and filtration efficiency.

- Exceptional inside filter for low breathing resistance
- Individually adjustable aluminum nose clip for excellent fit
- Ultrasonic welded head-straps for secure seal
- Optional one-way valve for easy exhalation
- High density urethane foam nose piece for great fitness and worker comfort

Standards

This product has been tested and meets the NIOSH standard of N95. Approved by NIOSH. The main performance tests in this standard are:

- Filter Penetration
- Breathing Resistance

N95 Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.

Uncle Mats and Supply, Inc.

Deerfield Beach, FL • 954-751-9800


info@unclesupply.com • www.UncleSupply.com



N-95 Particulate Respirators

Storage and Transportation

The N95 Series Particulate Respirators have a shelf life of 3 years. End of shelf life is marked on the product packaging. Before initial use, always check that the product is within the stated shelf life (use by date). Product should be stored in clean, dry conditions within the temperature range: -5°C to + 50°C with a maximum relative humidity of <75%. When storing or transporting this product use original packaging provided.

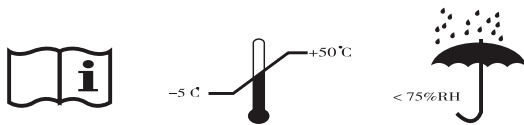

Limited Warranty

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Vanguard will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Vanguard products proved to be defectively manufactured. Except as provided above, Vanguard shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special or consequential, arising out of sale, use, misuse or the inability to use products by the user.

Service Time Limitation

- Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted. Do not use the same respirator again once removed.
- The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

Uncle Mats and Supply, Inc.

Deerfield Beach, FL • 954-751-9800

info@unclesupply.com • www.UncleSupply.com

© 2020 VGuard® is a registered trademark of Vanguard International Solutions, Inc.

Warnings

1. Wearer must read fitting instructions and limitations before use.
2. This respirator does not protect against the risk of contracting disease or infection.
3. Individuals highly sensitive to natural rubber latex may have an allergic reaction.
4. Wearer must do fit test either quantitatively or qualitatively by a qualified fit tester before using the respirator in a contaminated environment, as per OSHA guidelines, 29 CFR 1910.134.

RESPIRATOR FITTING INSTRUCTION

1



STEP 1:

Hold the respirator in hand with the nose piece at your fingertips, allowing the headbands to hang freely below your hand.

2



STEP 2:

Press the respirator firmly against your face with the nose piece on the bridge of your nose.

3



STEP 3:

Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.

4



STEP 4:

Using both hands, mold nose piece to the shape of your nose.

5



STEP 5:

To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension for better fit.

Please carefully follow these fitting instructions during each use to achieve proper fit.